

# CVIP3 – CRYOGENIC VIBRATION ISOLATOR PLATFORM 3



#### **Features**

- 3 DoF vibration attenuation inside the cryostat
- Placed between cold plate and experiment
- High insensitivity to cables and thermal braids
- High Q factor (no forced damping)
- Nonmagnetic
- Large central open aperture

### **Description / Applications**

The CVIP3 is part of the CVIP series, a general use passive vibration isolation platform that is placed directly on the cold plate and carries the experiment, thus reducing z, Rx and Ry vibrations into the experiment. For practicality and ease of use the cut-off frequency, above which vibrations are attenuated, is designed in the tens of Hertz range thereby accepting transmittance of the less critical low frequency vibrations but enabling a substantial payload and having reduced sensitivity for attached wiring.

Specifications	
General info	
Isolated axes	Z, Rx and Ry
Dimensions	See drawings below
Operational environmental conditions	20 mK to 375 K, ambient to UHV
Weight	1,3 kg
Frequencies and load capacity	
Payload	6 kg
Cut-off frequencies	Z: 42 Hz / Rx: 30 Hz / Ry 30 Hz (max payload, centered 35mm above platform)
Mechanical end stops	+/-0.5 mm
Materials	
Construction	Phosphor bronze












#### **Ordering Information** Available models CVIP<sub>3</sub> Cryogenic Vibration Isolator Platform 3 **Available Options** -HV Upgrade to High Vacuum compatibility -UHV Upgrade to Ultra High Vacuum compatibility Accessories None Mechanical and electrical information Contact Download 3D step files and manuals from: For quotations, specials, or engineering services, please contact us at: https://www.jpe-innovations.com/cryo-nano-products/ https://www.jpe-innovations.com/contact/

## Drawings

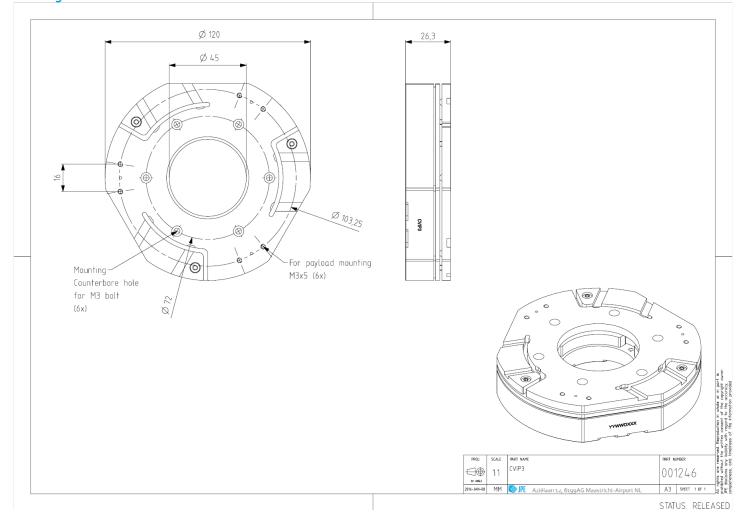

JPE



+31 (0) 43 3585 777



© JPE. Subject to change without notice. All data are outdated by any new release. JPE disclaims any liability with regard to the accuracy, completeness, and timeliness of the information provided. Brochure date: 17 February 2023, page 2 of 2